



# RISH SURGICAL

+92 3365240977

+92 3365240977

rishsurgicalinstruments@gmail.com

## Hysterectomy Set



## Rochester Ochsner Forcep 20 cm, STR



Rochester-Ochsner Forceps are hemostats that are helpful for clamping off blood vessels to avoid bleeding during procedures such as orthopedic surgery.

## Heaney Vaginal Clamp 19/21 cm CVD



Heaney Vaginal  
Clamp 19 cm CVD



Heaney Vaginal  
Clamp 21 cm CVD

Heaney Heaney Vaginal Clamp 19 cm CVD are ratcheted, finger ring forceps used for grasping uterine ligaments during vaginal hysterectomy procedures. They are available straight or curved with one or two teeth. Heaney forceps have diagonal serrations the full length of the jaws, and the jaws are approximately 1/3 the length of the shanks. The combination of serrations and teeth provide a secure hold on the ligaments. This product has a single tooth which is located slightly below the tip. It is curved and has an overall length of 195 mm. It also features a heavy pattern.

## Wertheim Vaginal Forcep 25 cm



Wertheim Vaginal Forcep is a forceps-like device that gynaecology and obstetrics specialists commonly use to manipulate heavier tissues, such as ligaments and other adnexa, during gynaecological surgery.

## Halsted Mosquito Forcep 12 cm, CVD/STR



Halsted Mosquito  
Forcep 12 cm, CVD



Halsted Mosquito  
Forcep 12 cm STR

Halsted-Mosquito Forceps is a specialized device that surgeons use to clamp blood vessels, especially small arteries, during a variety of surgical scenarios. The instrument is also useful to grasp disposable materials, tissues and sutures.

## **Crile Forcep 14 cm, 1:2 CVD**



Crile Forceps is a specialized device that surgeons use to clamp blood vessels, especially small arteries, during a variety of surgical scenarios. The instrument is also useful to grasp disposable materials, tissues and sutures.

## **Mayo Hegar Needle Holder 20 cm**



Mayo-Hegar Needle Holders are often used in general surgery for clamping onto suture needles of various lengths.

## **Doyen Intestinal Forceps 23 cm, CVD**




Doyen Intestinal Forceps is a multipurpose surgical instrument commonly used to grasp internal organs as well as secure bleeders during surgeries. It is often used to grip intestinal walls during complex general surgical procedures, as well as gynecological and cardiovascular surgeries.

## **Masson Needle Holder 27 cm**



Masson Needle Holders are used by surgeons to hold and mobilize needles during suturing procedures.

## **Allis Forceps 15 cm, 5x6**

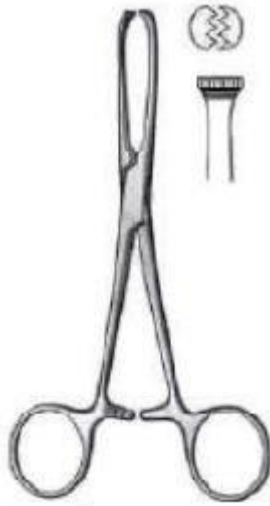

Allis Forceps is a popular surgical instrument that surgeons commonly use to grasp, statically hold or mobilize tissues of several kinds, in order to allow dissection, expose the surgical field or accessing a cavity.

## **Pozzi Tenaculum Forcep 25.5 cm**



Pozzi Tenaculum Forceps is a specialized OB/GYN device that surgeons can use for holding multiple structures, such as the cervix, uterus, blood vessels and the fallopian tubes, in order to perform multiple gynaecological procedures.



## **Metzenbaum Scissor 23 cm, BL/BL and Metzenbaum Scissor 25 cm, BL/BL, CVD**



**Metzenbaum Scissor  
23 cm, BL/BL**



**Metzenbaum Scissor  
25 cm, BL/BL, CVD**

Metzenbaum Dissecting Scissor is a versatile operating scissor ideal for dissecting delicate tissues in confined places. It is also used in general surgery for precise dissection.

## **Mayo Scissor 14 cm, BL/BL, STR and Mayo Scissor 17 cm, BL/BL, CVD**



**Mayo Scissor  
14 cm, BL/BL, STR**



**Mayo Scissor 17 cm,  
BL/BL, CVD**

Mayo Dissecting Scissor is a cutting device that surgeons in several fields use to sever superficial soft tissues, debris, muscles, fascia and deep heavy structures, in order to excise organs or to broaden the operating site.

## **Bozemann Forcep 26 cm, STR**



Bozeman Forceps is a specialized surgical tool that surgeons in multiple areas, especially in GYN/OB, can use to absorb blood and other fluids from the surgical field.

## **Schroeder Tenaculum Forcep 25 cm**



Schroeder Tenaculum Forceps are commonly used in surgery to hold blood vessels or steady the cervix and uterus. The forceps feature sharply pointed hooks at the tips of ring handles for comfortable and controlled use.



## **Ramsey Forcep 25 cm, S/J**



Ramsey Sponge Holding Forceps is a multipurpose surgical tool that surgeons frequently use to grab and hold swabs and sponges, which help in cleaning the blood and other fluids from the operating site.

## **Sims Scissor 20 cm, CVD**




Sims Gynaecological scissor is used for a variety of purposes such as cutting away bandages, trimming dead skin, dissect tissue and/or manipulate tissue if need be.

## **Spencerwells Forcep 15 cm CVD**



They are used during surgery to compress an artery, clamp and seal small blood vessels or to hold the artery out of the way.

## **Backhaus Towel Forcep 9 cm**

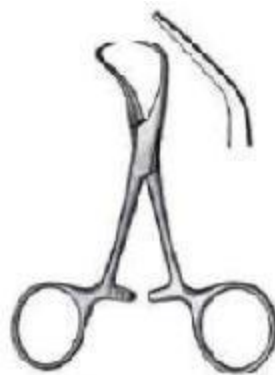


The Backhaus towel clamp is a perforating clamp. It is used for grasping the tissue, securing towels or drapes, and holding or reducing small bone fractures.

## **Doyen Myoma Instrument 17 cm**



Doyen Myoma Screw is a specialized gynecological instrument used to dig out and remove myomas from the uterine cavity.

## **Balfour Retractor 100 mm**

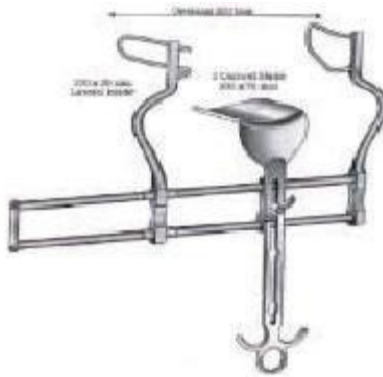

Balfour Abdominal Retractor is a self-retaining retractor used in laparotomy procedures. It may also be used for specific abdominal procedures where the abdomen needs to be held open for examination or evaluation, such as cesarean sections and bowel resection.

**Forceps Narrow Dress. 14.5 cm, Standard Forceps 1X2T 14.5 cm, Standard Dressing Forcep 25 cm, Standard Tissue Forceps 20 cm, 1:2**



Forceps Narrow Dress. 14.5 cm

Standard Forceps 1X2T 14.5 cm

Standard Dressing Forcep 25 cm

Standard Tissue Forceps 20 cm, 1:2

Forceps are used for holding dressing materials such as cotton and gauze during surgical procedures, changing dressings, or packing wounds

**Scalpel Handle No. 4**



Scalpel Handle No 4 is a hand-held device that surgeons in multiple fields use to carry and hold scalpel blades

## Scalpel Handle Standart Fig 3



Scalpel Handle No 3 is a specialized device that surgeons in several areas use to statically hold scalpel blades, in order to incise the skin with maximum precision

## Collin Duval Forceps 20 cm, 14 mm



Collin-Duval Intestinal Forceps is a specialized instrument often used for gripping bulky tissues. Its innovative tips minimize damage while ensuring a secure grip. Commonly used in general, gynecological, orthopedic procedures as well as other specialized surgeries.

## **Doyen Vaginal Specula 115x45 mm and Doyen Vaginal Specula 85x45 mm Fig 2**



Doyen Vaginal  
Specula 115x45 mm

Doyen Vaginal  
Specula 85x45 mm  
Fig 2

Doyen Vaginal Speculum is a specialized instrument used to visualize the upper genital tract during out patient procedures. It is often used during therapeutic and diagnostic gynecological procedures.

## **Doyen Retractor 28 cm, 48x90 mm and Doyen Retractor 28 cm, 48x105 mm**



Doyen Retractor  
28 cm, 48x90 mm

Doyen Retractor  
28 cm, 48x105 mm

Doyen Retractor is a multipurpose surgical tool that allows surgeons to widen surgical incisions for exposing their underlying cavities, hidden spaces, and the organs within the operating site.

## **Haberer General Instrument 30 cm 37/45**



Haberer Abdominal Spatula is a hand-held device that allows surgeons to pull back heavy tissues and the large organs inside the abdominal cavity in order to obtain maximum exposure of the surgical field.

## **Simon Retractor 28 cm 115x27 mm**



Simon Retractor is a hand-held surgical tool that allows surgeons to grasp, retract and maintain tissues away from the surgical field, especially during procedures involving the pelvic cavity or during gynecological surgery.



**Sterilization Container Stainless Steel 400x200x130 mm**

