



RISH SURGICAL

+92 3365240977

+92 3365240977

rishsurgicalinstruments@gmail.com

Mastoidectomy Set



Halsey Needle Holder Straight 130 mm (5 1/8")



Halsey Needle Holders are finger ring forceps that assist surgeons in guiding and inserting suture needles.

Delicate Scissors Straight 110 mm (4 1/4") SH/SH and Delicate Scissors Straight 120 mm (4 3/4") SH/SH



**Delicate Scissors
Straight 110 mm
(4 1/4") SH/SH**



**Delicate Scissors
Straight 120 mm
(4 3/4") SH/SH**

Delicate Operating Scissor is a multipurpose instrument that is used for sharp dissection of sensitive tissues. It is commonly used for dissecting organ tissue in Neuro, Plastic, ENT and General surgeries.

Mayo Stille Dissecting Scissors, Curved 140 mm (5 1/2") Chamfered Blades BL/BL



Mayo-Stille Dissecting Scissor is a sturdy general surgical instrument designed for deep dissection of dense/tough tissues during general surgeries.

Providence Hospital Hemostatic Forceps Curved 140 mm (5 1/2") Delicate, Blunt



Providence-Hospital Hemostatic Forceps is a scissor-like device that surgeons in different areas can use to manipulate blood vessels, clamp their lumen and control the blood flow, in order to achieve adequate haemostasis.

Ramsey Sponge Forceps Straight 180 mm (7") Serrated Fenestrated, Box Lock with Ratchet



Ramsey Sponge Forceps is a multipurpose surgical tool that surgeons frequently use to grab and hold swabs and sponges, which help in cleaning the blood and other fluids from the operating site.

Towel Clamp, Straight 85 mm (3 3/8") English Pattern, Sharp




Towel Clamp is a commonly used surgical instrument that is used to grasp tissues, drapes, gauze, and towels during surgeries.

Hartmann Nasal Speculum Med 160 mm



Hartmann Nasal Speculum assists doctors in endoscopic nasal procedures. They retract nasal walls and allow inspection of the cavities for diseases.

Tilley Ear Tampon Forceps Straight 140 mm (5 1/2") Knee Bent Very delicate

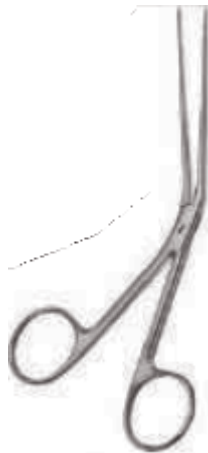

Tilley Ear Tampon Forceps are scissor handle forceps with serrated jaws that assist otolaryngologists to grasp and remove polyps in the ear.

**Volkman Wound Retractor (with Prong) 220 mm (8 3/4”),
4 Prongs, SH, Depth 9 mm, Width 19 mm**



Volkman Retractor is a popular hand-held instrument that allows surgeons to manipulate tissues during distal limb surgeries, especially those involving small joints, soles, and palms.

Jansen Bone Curette 160 mm (6 1/4”) Double Ended



Jansen Bone Curette is an orthopedic surgical device that surgeons commonly use for scraping fibrous, bony and diseased tissue from several bone structures, as well as soft tissue entrapped between fracture edges.

Scalpel Handle 160 mm (6 1/4") No. 7



Scalpel Handle No. 7 is designed to be used with scalpel blades No. 10 to 17

Weitlaner Retractor (Self Retaining) 165 mm (6 1/2") 3X4 Prongs, BL and Weitlaner Retractor (Self Retaining) 130 mm (5 1/8") 3X4 Prongs, BL



Weitlaner Retractor
(Self Retaining) 165 mm
(6 1/2") 3X4 Prongs, BL



Weitlaner Retractor
(Self Retaining) 130 mm
(5 1/8") 3X4 Prongs, BL

Weitlaner Self Retaining Retractor is a specialized instrument commonly used to retract incisions for better surgical access. It does so with minimal damage to the underlying tissues. Commonly used in general, neurological, spinal and orthopedic surgeries.

Troeltsch Micro Forceps Nerve Root Retractor Straight 115 mm (4 1/2") Kneebent Serrated Width 1.80 mm



Troeltsch Micro Forceps Nerve Root Retractor is a specialized device that neurosurgeons can use for pulling back, hold and mobilize spinal nerve roots and other soft tissues during neurosurgical procedures.

Thudichum Nasal Speculum #2



Thudichum Nasal Speculum assists doctors to examine and treat nasal diseases. It is routinely used in diagnostic and therapeutic nasal procedures, like anterior rhinoscopy for foreign body removal, nasal packing and surgeries of the nasal septum.

Gillies Tissue Forceps Straight 155 mm (6") Delicate Toothed (1X2)



Gillies Dressing Forceps is a multipurpose surgical tool that surgeons in every area can use for grasping swabs, drapes, gauzes and other small-sized surgical dressings, with the purpose of packing wounds and absorbing excess fluids

MC Indoe Dressing Forceps (Tweezers) Straight 150 mm (6") Delicate, Serrated



Mc Indoe Dressing Forceps is a multipurpose surgical device that surgeons in several specialities can use to grasp, hold and mobilized either tissues and dressings, in order to maintain a clear view of the surgical field.

Gouge 4, 6, 8, 10, 12 mm 140 mm 5 1/2"



Gouge is a specialized orthopedic device that surgeons can use to cut cortical bone tissue or to scoop the soft, cancellous bone collected for bone graftings.

Gillies Skin Hook 180 mm (7") Small 1 Prong SH



Gillies Skin Hook is an extremely useful, multipurpose device, which allows for grasping, holding, and elevating flaps of the skin and other soft tissues, especially during wound exploration, biopsy harvesting, and reconstructive surgery procedures.

Howarth Elevator Curved 215 mm (8 1/2") Double Ended SH/BL



Howarth Periosteal Elevator are double-ended instruments that are used in the elevation of the mucoperichondrium and the retraction of soft tissues in nasal surgeries.

Mallet 185 mm (7 1/4") Dia. 25 mm, Stainless Steel 210G, 135G



Bone Mallet is a very common orthopedic surgical instrument that surgeons can use for hammering the impacting platforms of bone chisels, tampers, and osteotomes, in order to perform osteotomy procedures and harvest bone grafting.

Hartmann Grasping Forceps (Small) Straight 80 mm (3 1/4") Serrated 1X8 mm



Hartmann Grasping Forceps is a multipurpose ENT tool that surgeons can use for grasping and dissecting tissues, packing the middle ear and holding small materials during otology procedures.

Sterilization Container SS

